

120 Drive with 400D1 Pumphead



100 series cased pump

Features and benefits

- Smooth, wipe down case with no sharp edges to tear gloves and no bug traps. IP31 ingress protection
- Superior speed control up to 2,000:1, $\pm 1\%$ set speed accuracy, repeatable flow
- Efficient and intuitive to use, requiring minimal key presses
- Compact footprint requires minimal countertop space
- Maintenance-free, brushless DC motor
- 1 channel pumphead
- Automatic variable speed, Manual variable speed options
- Three-year warranty



120 Drive with 400D1 Pumphead performance

400D1 pumphead for continuous tubing for manual (120S) and auto (120U) variable speed pumps — Tube bore and flow rates (ml/min)							
	Speed rpm	0.5mm	0.8mm	1.6mm	2.4mm	3.2mm	4.0mm
120S/D1	1 - 200rpm	0.01 – 2.2	0.03 – 5.8	0.11 – 23	0.24 – 49	0.41 – 81	0.59 – 120
120U/D1	0.1 - 200rpm	0.001 – 2.2	0.003 – 5.8	0.011 – 23	0.024 – 49	0.041 – 81	0.059 - 120

Technical specifications


	120 Drive with 400D1 Pumphead
Pumphead number of channels	1
FlowRateRange	0.001 ml/min to 120 ml/min
FlowRateRange	0.00002 USGPH to 3.04 USGPH
Max. operating pressure	30 psi (2 bar)
Speed control ratio	200:1, 2000:1
	0.1 rpm to 200
Operating temperature	40 °F to 104 °F (5 °C to 40 °C)
Weight	3.6 lbs (1.65 kg)
Control options	Automatic variable speed, Manual variable speed
Drive control options	S, U
Standards	EN 61010-1:2010/A1:2019, IEC 61010-1:2010/AMD1:2016
Drive ingress protection	IP31, NEMA 2
Drive humidity	(Non-condensing) 80% up to 31°C (88°F), decreasing linearly to 50% at 40C (104F)
Noise	<60dB(A) at 1 m
Drive noise	<60dB(A) at 1 m
Max. altitude	6562 ft (2000 m)
Drive power supply	Supplied with 24 VDC mains transformer plug (range 90 to 264 VAC 47 to 63 Hz) which comes with multi-country adaptors

Materials of construction

	120 Drive with 400D1 Pumphead
Driveshaft	Stainless steel
Drive casework	ABS plastic
Drive keypad/HMI	Polyester
Pumphead body assembly	Aluminum (anodized)
Pumphead roller assembly	Acetal
Pumphead rotor assembly	Aluminum (anodized)
Pumphead track	Aluminum (anodized)
Spindles	Stainless steel

Information listed covers the complete range.
For detailed specifications of individual models/components refer to user manual or contact WMFTS representative.

120 Drive with 400D1 Pumphead dimensions

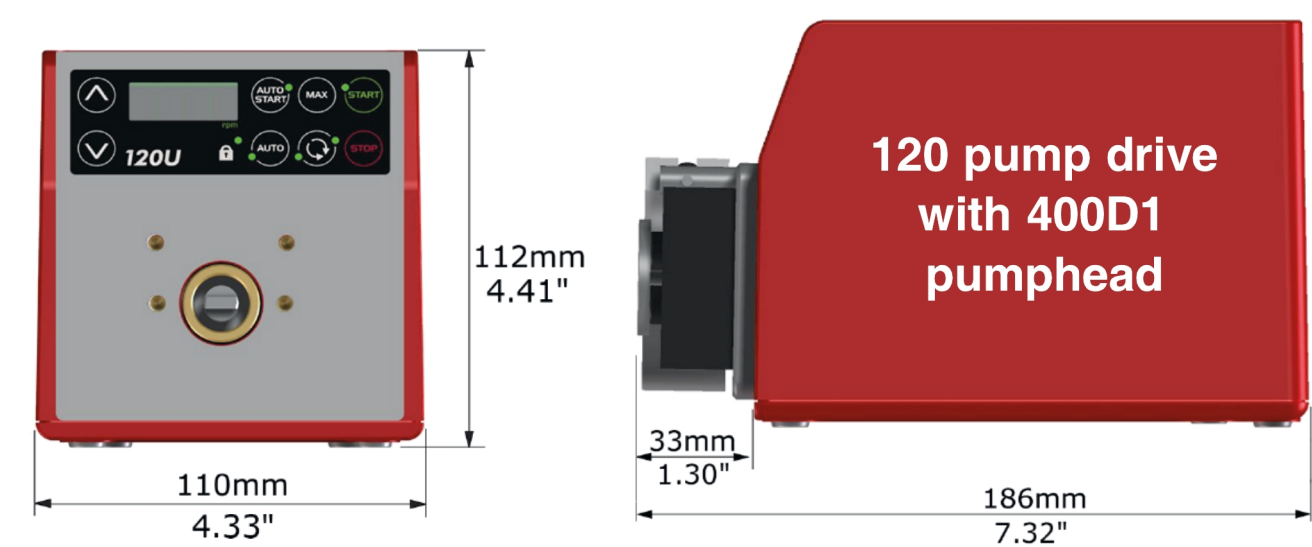


Control options

Key to pump drive control											
120F		120S					120U				
Fixed speed		Manual variable speed					Auto variable speed				
Manual control											
Pump drives							120F	120S	120U		
Fixed speed (CW direction)							✓				
Variable speed 1rpm increments								✓			
Variable speed 0.1rpm increments									✓		
Numerical speed display								✓	✓		
Direction change								✓	✓		
Max (prime)								✓	✓		
Auto start								✓	✓		
Remote control and feedback											
Pump drives							120F	120S	120U		
Run/stop, direction change, auto/man toggle, alarm									✓		
Analogue control											
Pump drives							120F	120S	120U		
Variable speed 0.1rpm increments									✓		
Analogue speed input; 4-20mA, 0-10V									✓		
Tacho speed output; 0-5V									✓		
Security											
Pump drives							120F	120S	120U		
Keypad lock to protect set-up								✓	✓		
Control Ratio											
	120F	120S					120U				
Pumphead	N/A	102R	114DV	400D1	400DM2	400DM3	102R	114DV	400D1	400DM2	400DM3
Manual	Fixed	1-32 rpm (32:1)	1-200 rpm (200:1)	1-200 rpm (200:1)	1-100 rpm (100:1)	1-100 rpm (100:1)	0.1-32 rpm (320:1)	0.1-200 rpm (2,000:1)	0.1-200 rpm (2,000:1)	0.1-100 rpm (1,000:1)	0.1-100 rpm (1,000:1)
Automatic	N/A	N/A	N/A	N/A	N/A	N/A	0.01-32 rpm (3,200:1)	0.01-200 rpm (20,000:1)	0.01-200 rpm (20,000:1)	0.01-100 rpm (10,000:1)	0.01-100 rpm (10,000:1)


Product codes

120F

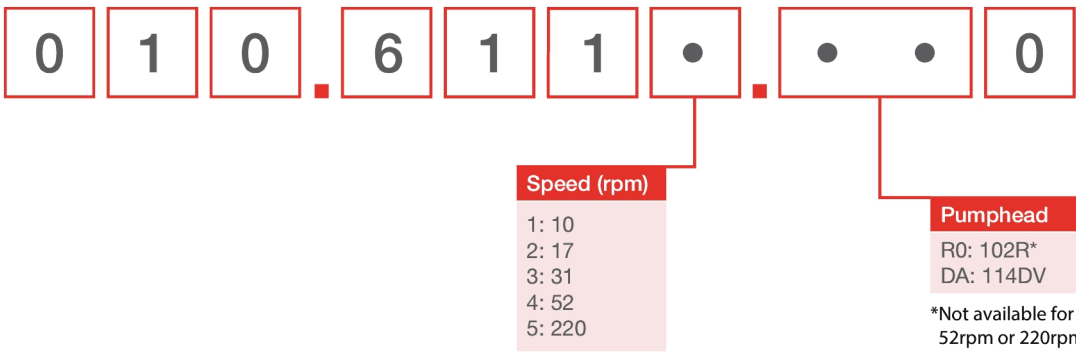

Note: All models are IP31 / NEMA 2

120S, 120U

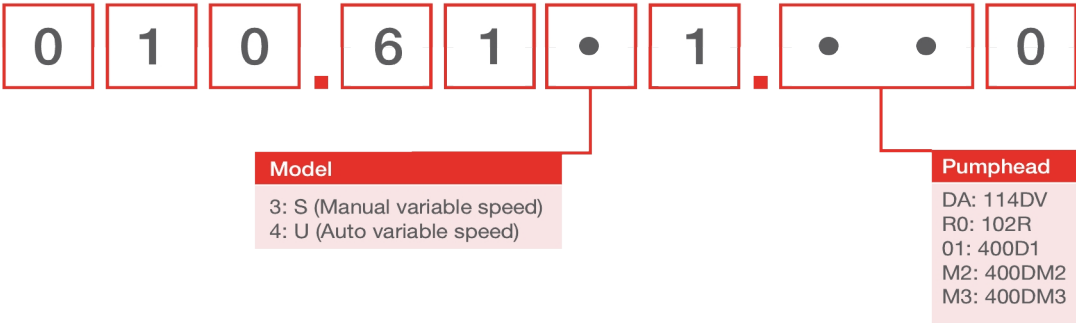

Pumphead product codes		
Model	Description	Product code
400D1	400D1 pumphead kit, 1 channel	043.001D.D1C
120 PSU	120V PSU kit 1A	MN2634B

Notes:
The pumpheads can only be sold as 120 replacement heads
The pumpheads are not interchangeable between drive models as they use different shaft designs

Tubing 1.6mm wall thickness for 400D1 pumphead						
	0.5mm	0.8mm	1.6mm	2.4mm	3.2mm	4.0mm
Bioprene	903.0005.016	903.0008.016	903.0016.016	903.0024.016	903.0032.016	
Pumpsil (Silicone)	913.A005.016	913.A008.016	913.A016.016	913.A024.016	913.A032.016	
STA-PURE® Series PCS			961.0016.016		961.0032.016	
Marprene	902.0005.016	902.0008.016	902.0016.016	902.0024.016	902.0032.016	902.0040.016
STA-PURE® Series PFL			966.0016.016		966.0032.016	
PVC			950.0016.016		950.0032.016	

Disclaimer: All flow rates shown were obtained pumping water at 20 °C (68 °F) with zero suction and delivery heads. The information contained in this document is believed to be correct but Watson-Marlow Limited accepts no liability for any errors it contains and reserves the right to alter specifications without notice. It is the users responsibility to ensure product suitability for use within their application. Watson-Marlow, LoadSure, Pumpsil, PureWeld XL, Bioprene, Marprene are registered trademarks of Watson-Marlow Limited. Tri-Clamp is a registered trademark of Alfa Laval Corporate AB. GORE and STA-PURE are registered trademarks of W.L. Gore and Associates. Please state the product code when ordering pumps and tubing.



wmfts.com/global
11 March 2026