

Maximise bioprocess efficiency



INNOVATION IN FULL FLOW

Market specific solutions for biotechnology and pharmaceutical applications

Watson-Marlow Fluid Technology Solutions delivers flexible and scalable fluid management solutions. Whether in single-use or stainless steel process, our range of Watson-Marlow peristaltic pumps, ASEPCO valves, BioPure single-use fluid path components, Aflex PTFE hoses, and Flexicon aseptic fill/finish systems all work together to reduce risk while optimising the repeatability of your process. All products offer traceability through LOT numbers or serial numbers.



Photo courtesy of NIBRT

WATSON MARLOW Pumps

Watson-Marlow peristaltic pumps and OEM solutions provide total process confidence and security, and scale from the benchtop to full production. Watson-Marlow's extensive product portfolio offers superior flow stability and metering accuracy ensuring process stability, GMP compliance and the highest final product quality.



100 series

Flow rates up to 190 ml/min. Single, twin or three channels of flow. Manual, remote or auto control.



200 series

Flow rates from 6 μ l/min to 22 ml/min. Cassettes accept standard manifold tubing. Up to 32 channels of flow.



300 series

Flow rates from 2 μ l/min to 2000 ml/min. Flip-top pumphead for easy tube loading. Single or multi-channel pumpheads.



500 series

Flow rates from 0.4 μ l/min to 3.5 L/min. Four drive and five pumphead options. Up to 7 bar pressure.



600 series

Flow rates from 0.001 ml/min to 18 L/min. Up to 4 bar pressure. 2650:1 control ratio.



700 series

Flow rates from 0.12 L/min to 55 L/min. Up to 2 bar pressure. Single or twin channels of flow.



800 series

Flow rates to 8000 litre/hr. 3.5 bar pressure. Clean-in-place (CIP) and steam-in-place (SIP).



Panel Mount Pumps

Flow rates from 0.01 μ l/min to 35L/min. Complete OEM pump with EtherNet/IP™ and PROFINET® options. Available with 100, 300, 400RXMD and 500 series pumpheads.



Flexicon is established as the preferred choice for aseptic liquid filling for GMP regulated industries including biotechnology and diagnostics. Our range of aseptic liquid filling equipment scales with your business, from stand-alone units for hand filling, through semi-automatic, to fully automatic filling, stoppering and capping machines.



Bench top fillers and cappers

Stand-alone units for manual aseptic filling and capping.



Semi-automatic systems

Semi-automatic highly flexible filling systems. Ideal for small batch production.



Automatic filling and capping

For clinical trial or contract filling. For bottles with screw caps or vials with crimp caps.



Fully automatic fill/finish system

The FPC60 is a uniquely configurable small batch aseptic fill/finish system.



OEM filling units

Peristaltic solutions for high capacity filling machines. Replacements for piston fillers.



Accusil® tubing and filling nozzles

Platinum-cured silicone tubing for accurate filling.



ASEPCO Weirless Radial diaphragm™ valves deliver consistency in biotechnology and pharmaceutical processes while reducing maintenance time by up to 80%. Their unique construction and radial diaphragm design makes them completely drainable and virtually eliminates the risk of contamination.



Weirless Radial diaphragm in-line valve family

Our family of Weirless Radial in-line valves are designed for a host of applications. In addition to the standard in-line valve, our sterile access valve and block-and-bleed valves are designed to eliminate the need for additional equipment and reduce cleaning and contamination risk. Our simple Tri-Clamp assembly allows for 80% faster maintenance.



Weirless Radial diaphragm tank bottom valve family

ASEPCO tank valves are also built around our patented radial diaphragm, and benefit from a simple Tri-Clamp assembly for faster maintenance. Tank valves can be welded truly flush to the bottom of your tank with no need for gaskets, bolts or seams. With our flexible design and manufacturing approach, we can customize valves to suit your needs.



Weirless Radial diaphragm sample valve

Our hydrodynamic design gives you clean, consistent sampling every time and allows for enhanced temperature control and sample system cleanability. The behind-the-seat flow path allows for easy CIP/SIP between samples. The design of the sample valve allows for closer jacketing to the tank with the fittings angled away from the vessel at a quicker rate.



BioPure produces a range of innovative bioprocessing fluid path components that simplify production operations, lower cGMP manufacturing costs and reduce process validation. All BioPure products are supported by class-leading validation guides detailing test protocols including ISO, USP and extractables, enabling customers to quickly assess suitability for their manufacturing process.



Q-Clamp™ sanitary Tri-Clamp™

Single-handed tool free operation. Unique tamper evident identification technology. Pressure rated up to 7 bar.



BioTube Applicator

Allows automation of tube to hose barb assembly. Increases quality and reduces risk of junction failure.



Platinum-cured silicone gaskets

Precision engineered to achieve a smooth bore, contamination free fluid path. Cleanroom manufactured and packed.



Single-use fluid path components

Tri-Clamps, flow control valves, end caps and connectors. Cleanroom manufactured and packed.

PolyClamp EPDM sanitary gaskets

SIP compatible with geometric stability after repeat cycles. USP Class VI validated.

PolyClamp EPDM gaskets are the perfect solution for connecting ASEPCO Weirless Radial diaphragm valves in your process.



WM architect single-use solutions

Optimize your process with a wide range of fluid management innovations from WMArchitect™ single-use solutions. Designed for biopharmaceutical and life science companies, WMArchitect will streamline your process, safeguard your product, reduce regulatory burden and resolve challenges for end-to-end, single-use fluid path management.



Customised designs

Optimize your process with custom-designed assemblies to meet your specific requirements. The open architecture model allows full flexibility and compatibility with existing fluid-contact materials used in your process.



Fluid Transfer

Optimize fluid transfer between equipment and process steps. Our specialist team has designed and standardized a selection of ready-to-use, pre-designed assemblies, all with lot traceable tubing and components.



Fill/Finish

Safeguard your therapeutic during the filling process with a wide range of sterility-assured assembly configurations. Representative assemblies are routinely tested for bioburden, particulates and endotoxin. Where further security is needed, our validation experts can provide lot-release testing to meet your quality requirements.

WATSON MARLOW Tubing

Watson-Marlow's range of high purity platinum-cured silicone, weldable and long life thermoplastic elastomer tubing is designed specifically for applications throughout biotechnology and pharmaceutical processes. Matching Watson-Marlow tubing with our peristaltic pumps enables customers to build easy-to-validate, end-to-end fluid path solutions from a single, reliable source.



Pumpsil®

Platinum-cured silicone with absolute traceability. Single-use tubing for biotechnology and pharmaceutical applications.



PureWeld XL®

Providing secure, cost effective pumping and transfer. Weldable, with low particulate generation.



Bioprene®

Thermoplastic elastomer tubing ideal for higher pressure applications with excellent chemical resistance.



GORE STA-PURE Series PCS

PTFE-reinforced silicone tubing with very low particulate generation. Delivers long term performance and repeatable accuracy.



GORE STA-PURE Series PFL

PTFE-reinforced fluoroelastomer tubing handles nearly all aggressive chemicals including organic solvents.




Platinum-cured silicone braided hose

High pressure flexible hose, single and double braided options. Cleanroom manufactured and packed.



Aflex hose manufacture the widest, most technically advanced range of PTFE lined flexible hose products. Performance features include excellent chemical resistance, a smooth bore that ensures a clean fast fluid flow, resistance to high pressures and vacuum resistant and unmatched flexibility and kink resistance.

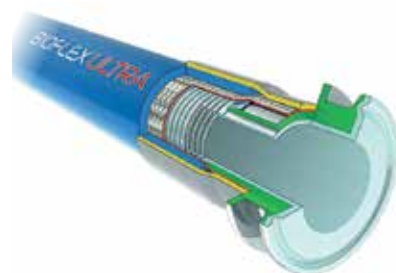

Bioflex Ultra

Design for process fluid transfer, Bioflex Ultra smooth bore hose has a patented liner that offers unmatched flexibility and kink resistance. Bioflex Ultra provides a clean, fast flow of high purity fluids. Bores available up to 3 inch (80mm) bore with a range of PTFE lined and non lined end fittings and design options.



Pharmaline N&X

Pharmaline N&X hose is designed to replace conventional silicone rubber hose for process fluid transfer. Pharmaline offers improved autoclave life cycle performance, chemical resistance and internal cleanability. Bores are available up to 3 inch (80mm) bore (Pharmaline N only) with a range of non lined end fittings and design options.



BIOTECHNOLOGY AND PHARMACEUTICAL SOLUTIONS



Watson-Marlow Fluid Technology Solutions

Watson-Marlow Fluid Technology Solutions supports its customers locally through an extensive global network of direct sales operations and distributors

[View full upstream range](#)



Disclaimer: The information in this document is believed to be correct, but Watson-Marlow Fluid Technology Solutions accepts no liability for any errors it contains, and reserves the right to alter specifications without notice. WARNING: These products are not designed for use in, and should not be used for, patient-connected applications. WMArchitect and Accusil are registered trademarks.