


ReNu SU Technology cartridge



Sterile ReNu SU Technology cartridge

Features and benefits


- Choose from a range of industry-standard aseptic connectors
- Maintain batch integrity
- Assemblies available with 10⁻⁶ sterility assurance level
- Maximize manufacturing flexibility and eliminate operator error
- Suitable for use with the Quantum pump
- Pressures up to 3 bar and flows up to 20 L/min

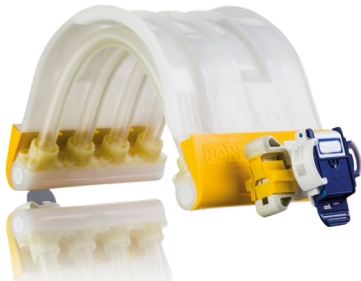


Technical specifications

	ReNu SU Technology cartridge
FlowRateRange	2 L/min to 20 L/min
FlowRateRange	120 L/h to 1200 L/h
Operating temperature	41 °F to 98.6 °F (5 °C to 37 °C)
Pump type	Cased
Connection	3/4" TC
Bore size	0.346 in (8.8 mm)
Shelf life	2 years
Static pressure limit	Do not exceed 6 Bar
Total cartridge volume	139 ml
Wetted surface area	Total area = 416 cm², Tubing (polyurethane) = 297 cm², Manifolds (HDPE USP VI) = 119 cm²
Priming	Self priming
Certification and compliance	ISO 10993-10, ISO 10993-11, ISO 10993-4, ISO 10993-6, USP <788>, USP <85>, USP <87>, USP <88>
BSE/TSE safe	Yes, compliant with EMEA 410/01 Rev 3
Tube gas Q permeability for oxygen	21 m²/(sPa)x10 ⁻¹⁸
Traceability	All material fully traceable back to raw materials
Dry running time	Do not exceed 2 mins
Weight	lbs (kg)


ReNu SU Technology cartridge dimensions

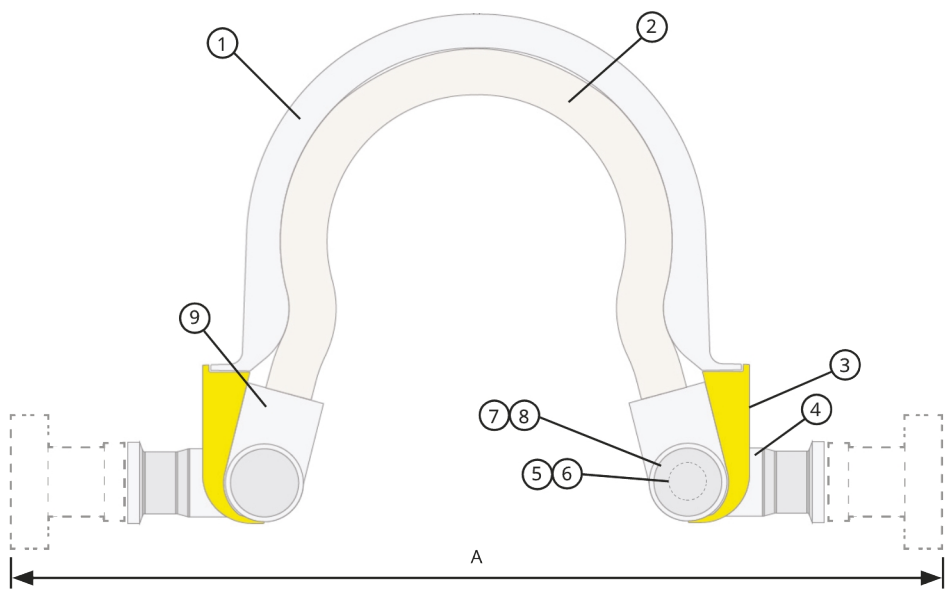


Standard Aseptic Connectors

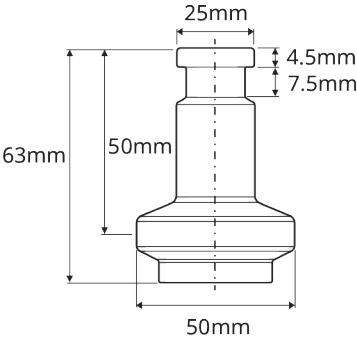

GE Readymate
(28-9366-95)

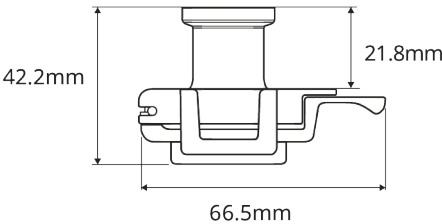

CPC AseptiQuik G
(AQG33012)

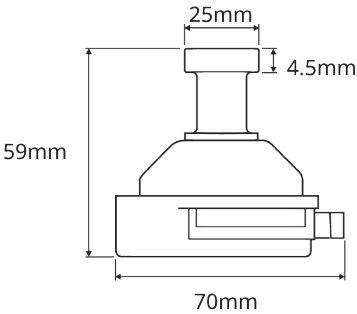

Pall Kleenpak Presto
(PSC1G08)

	Component	Material
1	Quantum track liner	HDPE (USPVI)
2	Quantum tube element	Polyurethane
3	Label, Product	Paper
4	Quantum manifold	HDPE (USPVI)
5	Quantum end cap	HDPE (USPVI)
6	Magnet	Neodymium, N35SH
7	Quantum cover plate	HDPE (USPVI)
8	Label, end cap	Paper
9	Quantum collet/joiner	Nylon 66
Aseptic connector		Dimension A (mm)
Kleenpak Presto Sterile Connector, 3/4" TC		313.0
Readymate Disposable Aseptic Connector, 3/4" TC		313.0
AseptiQuik G Genderless Connector, 3/4" TC		277.5
Cartridge only (no connectors)		195.0

Materials of construction - Cartridge assemblies and connection kits

Component	Material
Biopure Pt Silicone single braided hose	Platinum cured braid-reinforced silicone hose
BioBarb 5/8" Mini	Polypropylene
Q-Clamp Tamper Evident Cover, Yellow	Polypropylene
Q-Clamp, 1/2" - 3/4" TC	Glass Reinforced Nylon
Sanitary Gasket, 3/4" Unflanged	Silicone
BioEndCap	Polypropylene
Kleenpak Presto Sterile Connector, 3/4" TC	See manufacturers data - Pall (PSC1G08)
Readymate Disposable Aseptic Connector, 3/4" TC	See manufacturers data - GE (28-9366-95)
AseptiQuik G Genderless Connector, 3/4" TC	See manufacturers data - CPC (AQG33012)

Product codes

ReNu SU Cartridge Product Codes - RENU SU CARTRIDGE 600 20/3P

Partcode	Description	Additional details	End termination
33-1061-000001	RENU SU CARTRIDGE 600 20/3P	Not gamma irradiated	2 x 3/4" TC sealed with Bioendcap
33-1061-000002	RENU SU CARTRIDGE 600 20/3P, IRRADIATED	Gamma irradiated and released based on an exposure of 25 kGy to 45 kGy	
33-1062-000002	RENU SU 600 20/3P GAMMA WITH PRESTO IN/OUT-STR	Gamma irradiated at 25kGy to 45kGy and certified (SAL) level of 1 x 10-6 CFU	2 x Kleenpak Presto Sterile Connector, 3/4" TC (PSC1G08)
33-1062-000003	RENU SU 600 20/3P GAMMA WITH READYMATE IN/OUT-STR	Gamma irradiated at 25kGy to 45kGy and certified (SAL) level of 1 x 10-6 CFU	2 x Readymate Disposable Aseptic Connector, 3/4" TC (28-9366-95)
33-1062-000004	RENU SU 600 20/3P GAMMA WITH ASEPTIQUIK G IN/OUT-STR	Gamma irradiated at 25kGy to 45kGy and certified (SAL) level of 1 x 10-6 CFU	2 x AseptiQuik G Genderless Connector, 3/4" TC (AQG33012)

1. This assembly has been gamma irradiated and released based on an exposure of 25 kGy to 45 kGy.
2. A sterilisation validation study has been conducted according to the ISO11137-2, Method Vdmax25 to verify that that a master assembly will meet a sterility assurance level (SAL) of 10-6, when irradiated to a minimum of 25 kGy. This product has been released on the basis of being a section of the master assembly validation study.
3. Installation recommendations— Avoid sharp bends in pipework. Recommended minimum exit hose length is 2 meters.

STR Denotes that this product is supplied with sterility assurance

ReNu SU Technology Product Codes - Cartridge Assemblies

Partcode	Description	Assembly tubing length	Additional details	End termination
33-1063-000001	RENU SU 600 20/3P GAMMA ASEPTIQUIK G IN/OUT, 0.5/0.5 MTR STR	0.5m (x2)	Gamma irradiated at 25kGy to 45kGy and certified (SAL) level of 1 x 10-6 CFU	2 x AseptiQuik G Genderless Connector, 3/4" TC (AQG33012)
33-1063-000002	RENU SU 600 20/3P GAMMA ASEPTIQUIK G IN/OUT, 1.0/1.0 MTR STR	1.0m (x2)	Gamma irradiated at 25kGy to 45kGy and certified (SAL) level of 1 x 10-6 CFU	
33-1063-000003	RENU SU 600 20/3P GAMMA ASEPTIQUIK G IN/OUT, 2.0/2.0 MTR STR	2.0m (x2)	Gamma irradiated at 25kGy to 45kGy and certified (SAL) level of 1 x 10-6 CFU	
33-1063-000004	RENU SU 600 20/3P GAMMA ASEPTIQUIK G IN/OUT, 3.0/3.0 MTR STR	3.0m (x2)	Gamma irradiated at 25kGy to 45kGy and certified (SAL) level of 1 x 10-6 CFU	

ReNu SU Connection Kit Product Codes

Partcode	Description	Assembly tubing length	Additional details	End termination
33-1069-000001	RENU SU CONNECTION KIT, BRAID TUBING, 3/4" TC, 0.5 MTR-NON	0.5m	NA	2 x 3/4 TC sealed with Bioendcap
33-1069-000002	RENU SU CONNECTION KIT, BRAID TUBING, 3/4" TC, 1 MTR-NON	1.0m	NA	
33-1069-000003	RENU SU CONNECTION KIT, BRAID TUBING, 3/4" TC, 3 MTR-NON	3.0m	NA	
33-1069-000004	RENU SU CONNECTION KIT, BRAID TUBING, 3/4" TC, 2 MTR-NON	2.0m	NA	

Kits comprise of:

1 x tubing (Biopure 15.9mm ID Pt Silicone single braided hose) ends terminated with BioBarb™ 5/8" Mini.

Fittings (supplied loose):

2 x Q-Clamp Tamper Evident Cover, Yellow

2 x Q-Clamp, 1/2" - 3/4" TC

2 x Sanitary Gasket, 3/4" Unflanged

Packing/Sterility details: All components are manufactured and assembled in ISO 14644 class 7 (class 10,000, Grade C) cleanrooms of production facilities operating ISO 9001 quality management systems. The cleanrooms are air conditioned to ensure temperature and humidity.

Disclaimer: The information contained in this document is believed to be correct but Watson-Marlow Limited accepts no liability for any errors it contains and reserves the right to alter specifications without notice. It is the users responsibility to ensure product suitability for use within their application. Watson-Marlow, LoadSure, Qdos, ReNu, LaserTraceability, Pumpsil, PureWeld XL, Bioprene, Marprene are registered trademarks of Watson-Marlow Limited. Tri-Clamp is a registered trademark of Alfa Laval Corporate AB.

GORE and STA-PURE are registered trademarks of W.L. Gore and Associates.



wmfts.com/global
30 October 2025